



NatureBright[®]

CLARITY BULB REPLACEMENT INSTRUCTIONS

WARNING: RISK OF ELECTRIC SHOCK!


Unplug unit before servicing.

1. Place the Clarity™ Skincare Light Therapy on its face and remove the four (4) face-plate screws.

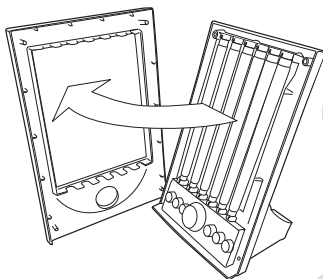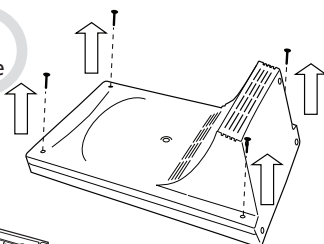


2. Remove cover carefully. Lay face-plate aside.



3. Turn the lamp so you can easily access the bulbs. Twist each bulb until the pins line up with the opening. Carefully remove bulbs by tipping the top of the bulb upwards. Then pull the bottom of the bulbs out of the fixture.

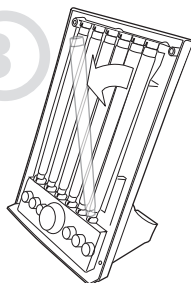

4. Insert a new bulb into the socket. Twist ¼ turn until the bulb locks in place.

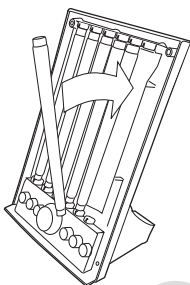

5. Replace cover.

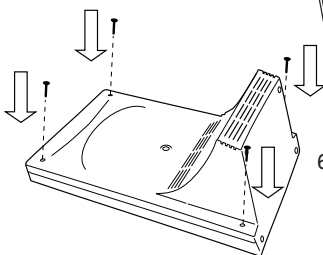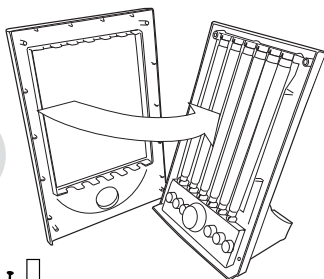

6. Lay the unit on face and replace four (4) screws. Tighten to snug fit.

Nature Bright Company 1-800-622-0231 www.naturebright.com
2600 Michelson Drive 17th Floor, Irvine, CA 92612
©2006 Nature Bright Company

Clarity_BulbReplace_0606